

FACE & BODY ROUTINES







FACE ROUTINES

INSTRUCTIONS FOR THE FACE	P3
RADIANCE ACTION	
MY RADIANCE ROUTINE	P4
ANTI-AGING ACTION	
MY EYES SMOOTHING ROUTINE	P5
MY ANTI-AGING FOREHEAD ROUTINE	P6
MY LIP CONTOUR ROUTINE	
MY ANTI-AGING NECK ROUTINE	P8
MY ANTI-AGING DÉCOLLETÉ ROUTINE	P9
MY ANTI-AGING HAND ROUTINE	P10
MY ANTI-AGING FACE ROUTINE	P11
FIRMNESS ACTION	
MY FACE FIRMING ROUTINE	P12
DETOX ACTION	
MY REFRESHED EYES ROUTINE	P13
MY DETOX ROUTINE	P14
SLIMMING ACTION	
MY FACE RESCULPTING ROUTINE	P15



Please carefully read the manual before using your Wellbox® [S].



> INSTRUCTIONS FOR THE FACE



The suction power must be adjusted depending on how you feel and on your skin quality. Be careful not to increase it too much on thin and mature skin.



These two pictograms are provided for information purposes only. They will guide you through the routine and indicate **which treatment head to use** according to the suction sequentiality you have selected.



Suction sequentiality is the number of aspirations per second. The higher the suction sequentiality, the closer to the skin surface will be the action; firmess or radiance.

Likewise, lower sequential suction will have a more in-depth action, ideal to fill in wrinkles.



For better tissue grab, **move the LIFT head following a dotted line** (with small hops, stopping 2-3 seconds on each point) while respecting the direction indicated by the white arrow.



The red lines indicate how to **position the flaps** when moving the treatment head. They must always be perpendicular to the wrinkle.



endermologie® cosmetics help to potentiate Wellbox®[S]'s action. Once stimulated, the skin becomes more receptive to the active ingredients contained in the cosmetics.



RADIANCE ACTION

> MY RADIANCE ROUTINE






STEP 1: I clean and prepare my skin with the **Pre-Care Micellar Lotion**, which soothes and provides maximum moisturization to my skin.



STEP 2: RADIANCE CARE





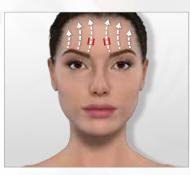












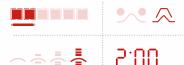



STEP 3: The **Total Eye Care** smoothes puffiness, wrinkles and dark circles from the eye **contour area**, my eyes are immediately refreshed.



STEP 4: The **Radiance Brightening Serum**, enriched with Vitamin C and the exclusive **LPG® Anti-Aging Complex**, brightens and evens my skin. It also can diminish my pigmentation spots.



STEP 5: Finally I apply the **Complete Correction Fluid SPF15** for an instantly-healthy glow. Lightly tinted, this fluid can smooth my skin's imperfections. Its city sun block (SPF 15) prevents the signs of aging by protecting from daily stresses.







> MY REFRESHED EYES ROUTINE 6 !!!

STEP 1: The **Pre-Care Micellar Lotion** should be used as a make-up remover/cleanser before each Wellbox[®][S] treatment. The gorse water soothes my skin and brings me maximum hydration.



STEP 2: REFRESHED EYES ROUTINE

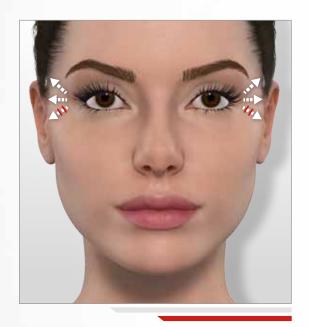

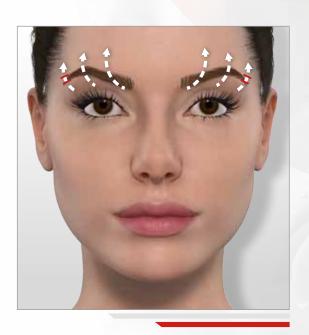

() [1]









~ * * * *

3:00

STEP 3: The **Total Eye Care** smoothes puffiness, wrinkles and dark circles from the **eye contour area**, my eyes are immediately refreshed.






> MY ANTI-AGING FOREHEAD ROUTINE




STEP 1: I clean and prepare my skin with the **Pre-Care Micellar Lotion**, which soothes and provides maximum moisturization to my skin.




STEP 2: ANTI-AGING FOREHEAD CARE

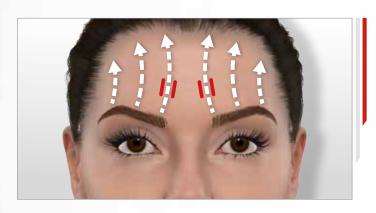



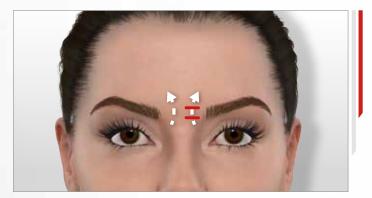

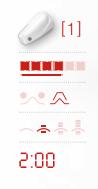

STEP 3: The **Radiance Brightening Serum,** enriched with Vitamin C and the exclusive LPG® Anti-Aging Complex, brightens and evens my skin. It also can diminish my pigmentation spots.



STEP 4: Thanks to its exclusive hyaluronic acid-based LPG® Anti-aging Complex, the **Replenishing Firming Cream** restores volume, smoothes and redensifies my skin.






> MY LIP CONTOUR ROUTINE




STEP 1: The **Pre-Care Micellar Lotion** should be used as a make-up remover/cleanser before each Wellbox[®][S] treatment. The gorse water soothes my skin and brings me maximum hydration.




STEP 2: LIP CONTOUR CARE

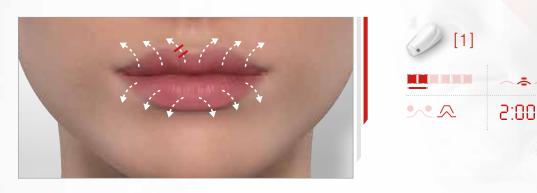

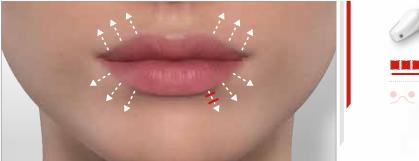



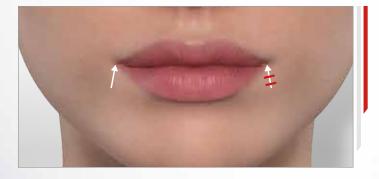



STEP 3: Thanks to its exclusive hyaluronic acid-based LPG® Anti-aging Complex, the **Replenishing Firming Cream** restores volume, smoothes and redensifies my skin.







> MY ANTI-AGING NECK ROUTINE

STEP 1: I clean and prepare my skin with the **Pre-Care Micellar Lotion**, which soothes and provides maximum moisturization to my skin.



STEP 2: ANTI-AGING NECK CARE

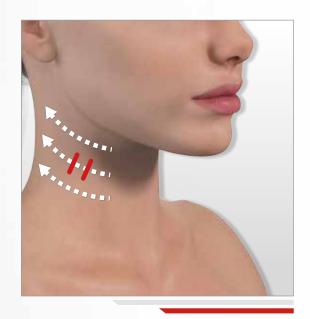

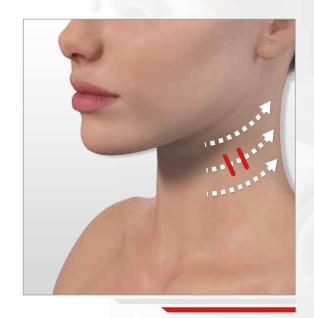

6 min













3:00

~ ~ **.** .

3:00

~ ~ \$ ****

STEP 3: The **Radiance Brightening Serum,** enriched with Vitamin C and the exclusive LPG[®] Anti-Aging Complex, brightens and evens my skin. It also can diminish my pigmentation spots.



STEP 4: Thanks to its exclusive hyaluronic acid-based LPG® Anti-aging Complex, the **Replenishing Firming Cream** restores volume, smoothes and redensifies my skin.







> MY ANTI-AGING DÉCOLLETÉ ROUTINE 6!


STEP 1: I clean and prepare my skin with the **Pre-Care Micellar Lotion**, which soothes and provides maximum moisturization to my skin.




STEP 2: ANTI-AGING DÉCOLLETÉ CARE

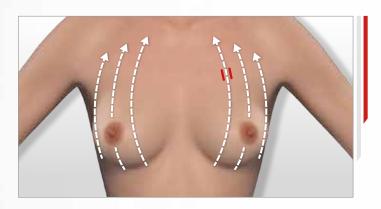

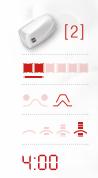

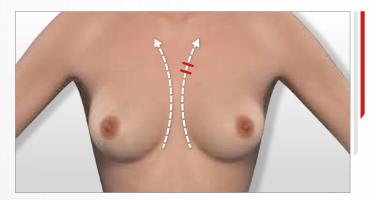

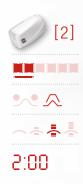

STEP 3: The **Radiance Brightening Serum,** enriched with Vitamin C and the exclusive LPG[®] Anti-Aging Complex, brightens and evens my skin. It also can diminish my pigmentation spots.



STEP 4: Thanks to its exclusive hyaluronic acid-based LPG® Anti-aging Complex, the **Replenishing Firming Cream** restores volume, smoothes and redensifies my skin.






> MY ANTI-AGING HAND ROUTINE 6 !!!

STEP 1: I clean and prepare my skin with the **Pre-Care Micellar Lotion**, which soothes and provides maximum moisturization to my skin.




STEP 2: ANTI-AGING HAND CARE



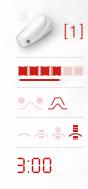



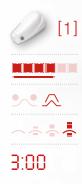

STEP 3: The **Radiance Brightening Serum**, enriched with Vitamin C and the exclusive LPG® Anti-Aging Complex, brightens and evens my skin. It also can diminish my pigmentation spots.



STEP 4: Thanks to its exclusive hyaluronic acid-based LPG® Anti-aging Complex, the **Replenishing Firming Cream** restores volume, smoothes and redensifies my skin.





> MY ANTI-AGING FACE ROUTINE






PREMIUM

STEP 1: The **Pre-Care Micellar Lotion** should be used as a make-up remover/cleanser before each Wellbox[®][S] treatment. The gorse water soothes my skin and brings me maximum hydration.



STEP 2: ANTI-AGING FACE CARE

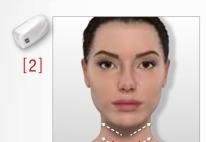



■■■■■■ 3:00



2:00





2:00

The Total Eye Care smoothes puffiness, wrinkles and dark circles from the eye contour area, my eyes are immediately refreshed.



STEP 4: The **Anti-aging Renewal Serum** concentrates 22% of active ingredients recognized for their anti-aging power. They act at the heart of the fibroblast to increase the production of collagen and elastin. This serum refines my skin, makes it brighter, better hydrated and my wrinkles are smoothed.



STEP 5: Thanks to its exclusive hyaluronic acid-based LPG® Anti-aging Complex, the **Replenishing Firming Cream** restores volume, smoothes and redensifies my skin.





FIRMNESS ACTION


> MY FACE FIRMNESS ROUTINE




STEP 1: I clean and prepare my skin with the **Pre-Care Micellar Lotion**, which soothes and provides maximum moisturization to my skin.



STEP 2: FACE FIRMNESS CARE

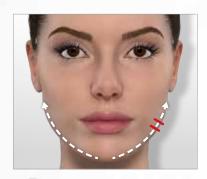


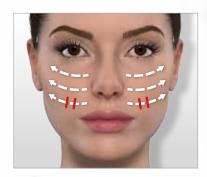







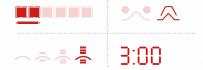

STEP 3: The **Total Eye Care** smoothes puffiness, wrinkles and dark circles from the eye contour area, my eyes are immediately refreshed.



STEP 4: The **Radiance Brightening Serum,** enriched with Vitamin C and the exclusive LPG® Anti-Aging Complex, brightens and evens my skin. It also can diminish my pigmentation spots.



STEP 5: Thanks to its exclusive hyaluronic acid-based LPG® Anti-aging Complex, the **Tensing Firming Cream** provides an immediate lifting effect. It also visibly reshapes facial contours.






DETOX ACTION


> MY REFRESHED EYES ROUTINE 6 min


STEP 1: The **Pre-Care Micellar Lotion** should be used as a make-up remover/cleanser before each Wellbox[®][S] treatment. The gorse water soothes my skin and brings me maximum hydration.




STEP 2: REFRESHED EYES CARE

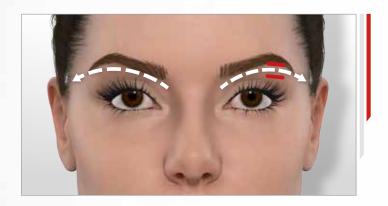

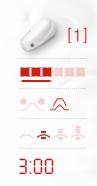

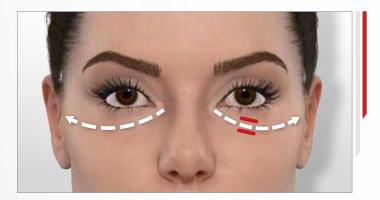

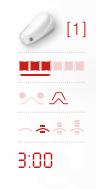

STEP 3: The **Total Eye Care** smoothes puffiness, wrinkles and dark circles from the eye contour area, my eyes are immediately refreshed.





DETOX ACTION

> MY DETOX ROUTINE






STEP 1: I clean and prepare my skin with the **Pre-Care Micellar Lotion**, which soothes and provides maximum moisturization to my skin.



STEP 2: DETOX CARE



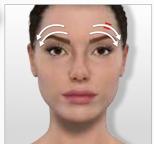

[1]



2:00 2:00





○♣♣♣ 2:00





~**♣**♣₩ 2:00



The **Total Eye Care** smoothes puffiness, wrinkles and dark circles from the eye contour area, my eyes are immediately refreshed.



STEP 4: Thanks to the **Intense Hydrating Smoothing Serum,** my skin is durably moisturized and regains its luscious and smooth quality. Its hyaluronic acid film also protects my skin from external stresses.



STEP 5: Finally I apply the **Complete Correction Fluid SPF15** for an instantly-healthy glow. Lightly tinted, this fluid can smooth my skin's imperfections. Its city sun block (SPF 15) prevents the signs of aging by protecting from daily stresses.






> MY FACE RESCULPTING ROUTINE 6 min

STEP 1: I clean and prepare my skin with the **Pre-Care Micellar Lotion**, which soothes and provides maximum moisturization to my skin.




STEP 2: FACE RESCULPTING CARE



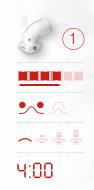



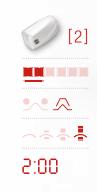

STEP 3: The **Firming V-Shaping Cream** helps to release my double chin fat and allows me to firm my facial contours, thanks to the exclusive LPG® Anti-aging Complex.





BODY

ROUTINES

INSTRUCTIONS FOR THE BODY	P17
FIRMING ACTION	8
MY ARM FIRMING ROUTINE	P18
MY STOMACH FIRMING ROUTINE	P19
MY BUTTOCKS FIRMING ROUTINE	P20
MY THIGH FIRMING ROUTINE	P21
SLIMMING ACTION	
MY TONED ARMS ROUTINE	P22
MY FLAT STOMACH ROUTINE	P23
MY ANTI-LOVE HANDLES ROUTINE	P24
MY SLIMMER THIGHS ROUTINE	
MY SLIMMER KNEES ROUTINE	
MY SLIMMER CALVES ROUTINE	. P27
ANTI-CELLULITE ACTION	
MY ANTI-CELLULITE BUTTOCKS ROUTINE	P28
MY ANTI-CELLULITE THIGHS ROUTINE	P29
MY ANTI-SADDLEBAGS ROUTINE	P30
WELLBEING ACTION	
MY LIGHT LEGS ROUTINE	P31



Please carefully read the manual before using your Wellbox® [S].



> INSTRUCTIONS FOR THE **BODY**



The suction power must be adjusted depending on how you feel and on your skin quality. Be careful not to increase it too much on sensitive areas, such as your inner arms and inner thighs.


These two pictograms are provided for information purposes only. They will guide you through the routine and indicate **which treatment head to use** according to the suction sequentiality you have selected.

Suction sequentiality is the number of aspirations per second. The higher the suction, the closer to the skin surface will be the action; firmess or radiance.

Likewise, lower sequential suction will have a more in-depth action, ideal to release resistant fat.



The white arrows indicate the direction that the treatment head must follow. **When using the ROLL3 head**, use the trigger to reverse the direction of the rollers.

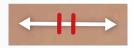

The red lines indicate how to **position the rollers** when moving the treatment head.



endermologie® cosmetics help to potentiate Wellbox®[S]'s action. Once stimulated, the skin becomes more receptive to the active ingredients contained in the cosmetics.




> MY ARM FIRMING ROUTINE 6 !!!

STEP 1: Every morning, I apply my **Glowing Resurfacing Body Cream,** for a new skin effect day after day. Its micro-peeling action acts effectively on my small skin imperfections (stretch marks, scars ...).



STEP 2: ARM FIRMING CARE

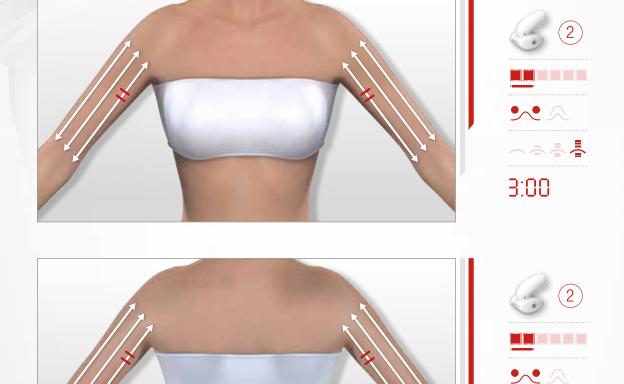

STEP 3: I apply the **Firming Contour Fluid** on the treated area for an immediate and lasting lifting effect.



3:00




> MY STOMACH FIRMING ROUTINE 6 11

STEP 1: Every morning, I apply my **Glowing Resurfacing Body Cream,** for a new skin effect day after day. Its micro-peeling action acts effectively on my small skin imperfections (stretch marks, scars ...).



STEP 2: STOMACH FIRMING CARE

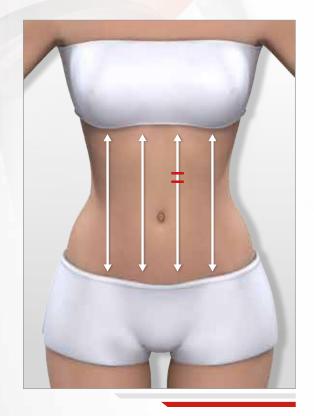















3:00

STEP 3: I apply the **Firming Contour Fluid** on the treated area for an immediate and lasting lifting effect.







> MY BUTTOCKS FIRMING ROUTINE 6 !!!


STEP 1: Every morning, I apply my **Glowing Resurfacing Body Cream,** for a new skin effect day after day. Its micro-peeling action acts effectively on my small skin imperfections (stretch marks, scars ...).




STEP 2: BUTTOCKS FIRMING CARE

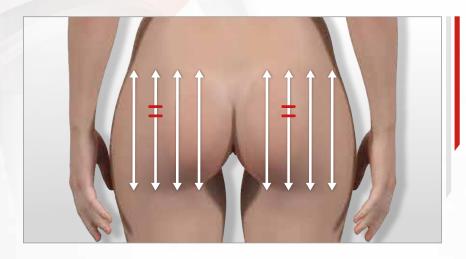

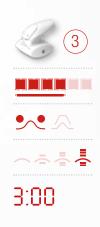

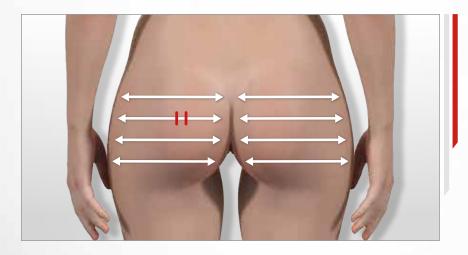

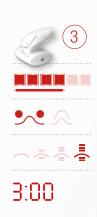

STEP 3: I apply the **Firming Contour Fluid** on the treated area for an immediate and lasting lifting effect.







> MY THIGH FIRMING ROUTINE 6 !!!

STEP 1: Every morning, I apply my **Glowing Resurfacing Body Cream,** for a new skin effect day after day. Its micro-peeling action acts effectively on my small skin imperfections (stretch marks, scars ...).



STEP 2: THIGH FIRMING CARE

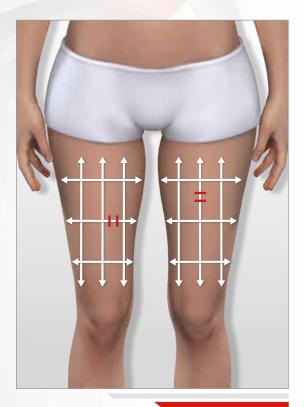

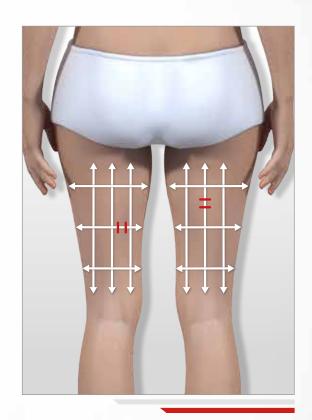













3:00

STEP 3: I apply the **Firming Contour Fluid** on the treated area for an immediate and lasting lifting effect.






> MY TONED ARMS ROUTINE 6 MI


STEP 1: Every morning, I apply my **Glowing Resurfacing Body Cream,** for a new skin effect day after day. Its micro-peeling action acts effectively on my small skin imperfections (stretch marks, scars ...).



STEP 2: TONED ARM ROUTINE

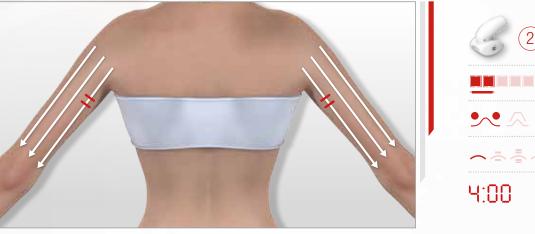



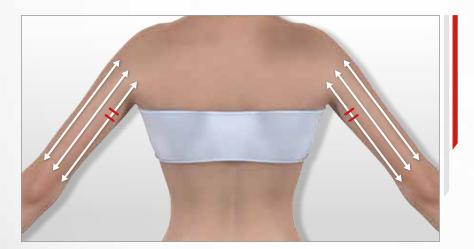







> MY FLAT STOMACH ROUTINE







STEP 1: Every morning, I apply my **Glowing Resurfacing Body Cream,** for a new skin effect day after day. Its micro-peeling action acts effectively on my small skin imperfections (stretch marks, scars ...).



STEP 2: FLAT STOMACH CARE

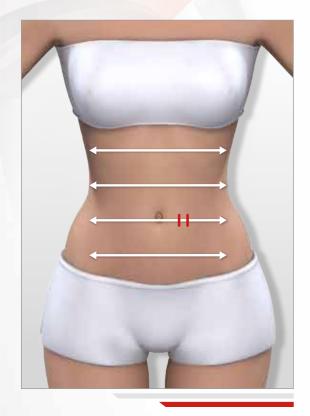

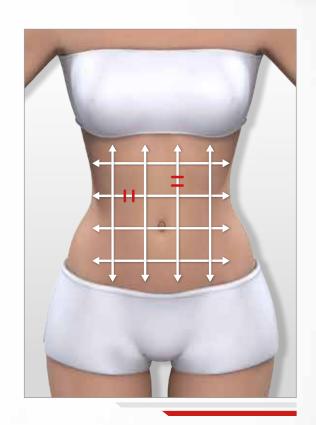













00:5

STEP 3: The exclusive LPG® Slimming Complex - containing pure caffeine, allows the **Body Shaping Gel** to reduce excess fat.





> MY ANTI-LOVE HANDLES ROUTINE



STEP 1: Every morning, I apply my **Glowing Resurfacing Body Cream,** for a new skin effect day after day. Its micro-peeling action acts effectively on my small skin imperfections (stretch marks, scars ...).



STEP 2: ANTI-LOVE HANDLES CARE



















STEP 3: The exclusive LPG® Slimming Complex - containing pure caffeine, allows the **Body Shaping Gel** to reduce excess fat.





> MY SLIMMER THIGHS ROUTINE 6!


STEP 1: Every morning, I apply my **Glowing Resurfacing Body Cream,** for a new skin effect day after day. Its micro-peeling action acts effectively on my small skin imperfections (stretch marks, scars ...).




STEP 2: SLIMMER THIGHS CARE



Very often slackening skin on the inside of the inner thighs is mistaken for a cellulite problem, when it actually is a lack of firmness. It is therefore not recommended to release fat or smooth this part of the body. You should instead you should opt for the **Firming Care.**

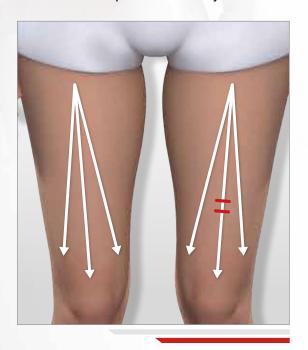

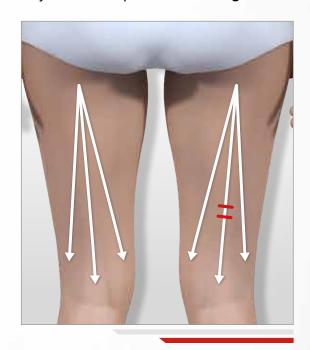











~ ~ * * *

3:00

It is highly recommended to firm the tissue after releasing fat. Think about finishing off your session with the THIGH FIRMING CARE.

STEP 3: The exclusive LPG® Slimming Complex - containing pure caffeine, allows the **Body Shaping Gel** to reduce excess fat.







> MY SLIMMER KNEES ROUTINE 6 !!!

STEP 1: Every morning, I apply my **Glowing Resurfacing Body Cream,** for a new skin effect day after day. Its micro-peeling action acts effectively on my small skin imperfections (stretch marks, scars ...).



STEP 2: SLIMMER KNEES CARE

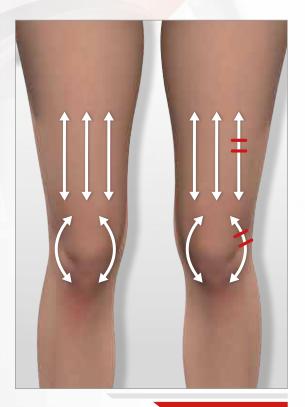

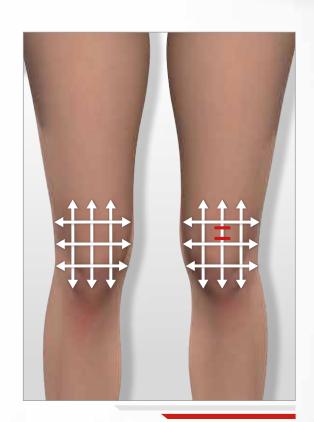











9:00

STEP 3: The exclusive LPG® Slimming complex - containing pure caffeine, allows the **Body Shaping Gel** to reduce excess fat.







> MY SLIMMER CALVES ROUTINE 6 min

STEP 1: Every morning, I apply my **Glowing Resurfacing Body Cream,** for a new skin effect day after day. Its micro-peeling action acts effectively on my small skin imperfections (stretch marks, scars ...).



STEP 2: SLIMMER CALVES CARE

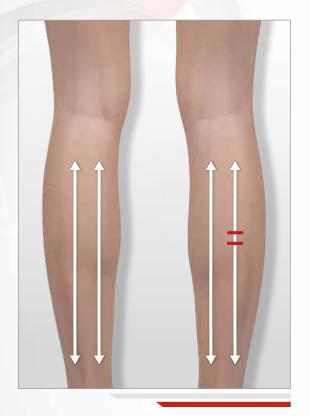

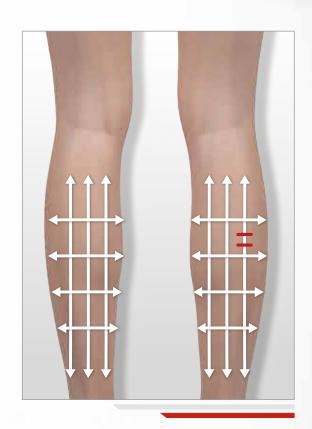











소축률 2:00

STEP 3: Enriched with the exclusive LPG® Slimming Complex, the **Draining Detox Serum** helps to quickly and effectively eliminate the roundness of the whole body. It is associated with Escin whose draining action effectively fights my aqueous cellulite.





ANTI-CELLULITE ACTION


> MY ANTI-CELLULITE BUTTOCKS ROUTINE





STEP 1: Every morning, I apply my **Glowing Resurfacing Body Cream,** for a new skin effect day after day. Its micro-peeling action acts effectively on my small skin imperfections (stretch marks, scars ...).




STEP 2: ANTI-CELLULITE BUTTOCKS CARE

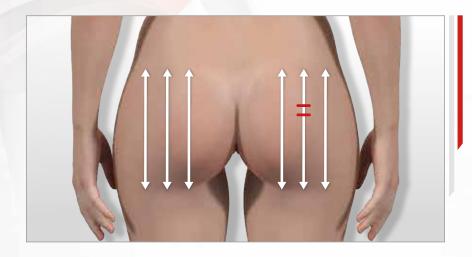

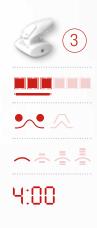

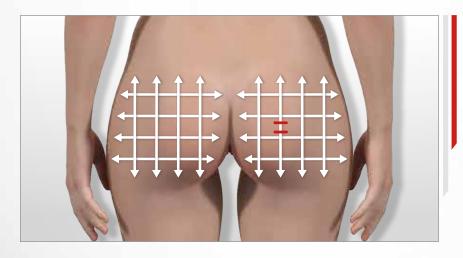

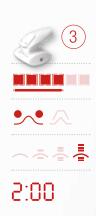

STEP 3: The exclusive LPG® Slimming Complex - containing pure caffeine, allows the **Orange Peel Smoothing Complex** to promote the elimination of fat and to smooth my cellulite.





ANTI-CELLULITE ACTION



> MY ANTI-CELLULITE THIGHS ROUTINE




STEP 1: Every morning, I apply my Glowing Resurfacing Body Cream, for a new skin effect day after day. Its micro-peeling action acts effectively on my small skin imperfections (stretch marks, scars ...).




STEP 2: ANTI-CELLULITE THIGHS CARE



Very often sagging skin on the inner thighs is mistaken for a cellulite problem, when it actually is a lack of firmness. It is therefore not recommended to release fat or smooth this part of the body. You should instead opt for the **Firming Care.**

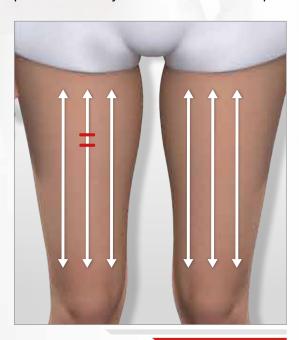

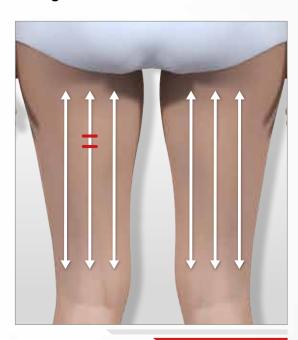











~~*

3:00

It is highly recommended to firm the tissue after releasing fat. Think about finishing off your session with the THIGH FIRMING CARE.

STEP 3: The exclusive LPG® Slimming Complex - containing pure caffeine, allows the **Orange Peel Smoothing Complex** to promote the elimination of fat and to smooth my cellulite.





ANTI-CELLULITE ACTION

> MY ANTI-SADDLEBAGS ROUTINE 6 miles

STEP 1: Every morning, I apply my **Glowing Resurfacing Body Cream,** for a new skin effect day after day. Its micro-peeling action acts effectively on my small skin imperfections (stretch marks, scars ...).



STEP 2: ANTI-SADDLEBAGS CARE




















2:00 ON EACH SIDE ~ ~ ~ **~**

CON EACH SIDE

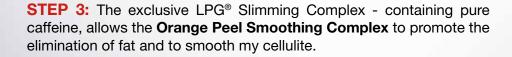





WELL-BEING ACTION

> MY LIGHT LEGS ROUTINE








STEP 1: Every morning, I apply my **Glowing Resurfacing Body Cream,** for a new skin effect day after day. Its micro-peeling action acts effectively on my small skin imperfections (stretch marks, scars ...).





STEP 2: LIGHT LEGS CARE

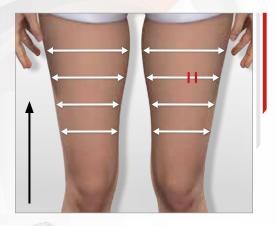

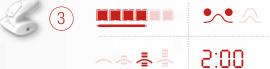

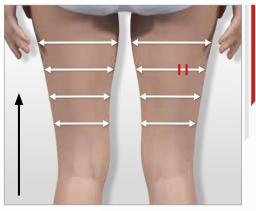

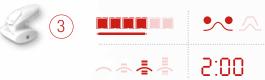

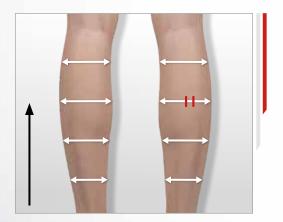



STEP 3: For maximum comfort, I spray the **Light Legs Mist** from the ankles upwards. It promotes stimulates micro-circulation and decongests my legs.





Find Wellbox® routines and all our scientific evidence on www.wellbox.com